

acid injection system



TEO LINK  
Syringe holder single use

TEOSYAL®PEN for your patients

- Significantly less pain
- Precise injections for a better result
- Reduced side effects

TEOSYAL®PEN manufactured by JuvaPlus



LABORATOIRES  
TEOXANE  
GENEVA

cdm  
paris

LET  
YOUR ART  
SPEAK



PRACTICAL  
GUIDE

TEOSYALPEN  
FOR THE FINE ART OF PRECISION

LABORATOIRES  
TEOXANE  
GENEVA

TEOSYAL®PEN. ART NOUVEAU  
The first motorized cordless hyaluronic

TEOSYALPEN

- Activation button
- On/Off
- Injection «Flow mode»  
3 Speeds: Slow/Medium/Fast
- Injection mode «Drop»  
3 Drop Sizes: Small/Medium/Large
- Rewind button




TEOSYAL®PEN for you


- Optimized control of the procedure:  
depth of injection, homogeneous delivery and precise dosage
- User-friendly & ergonomic, designed for more comfort
- Easy to use & intuitive device



1 Syringe preparation

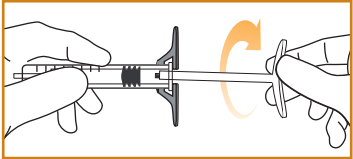


Unscrew the plunger rod from the syringe

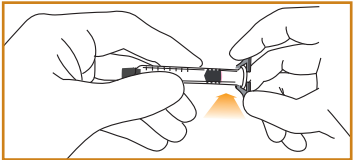

Remove and discard the finger grip from the syringe

TEOSYAL®PEN is compatible with the large majority of the hyaluronic acid brands


2 Setting up of the syringe

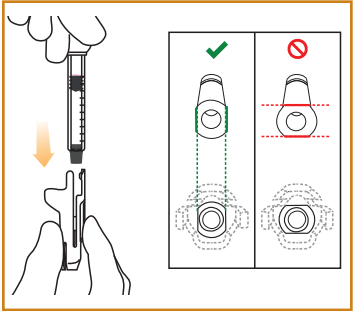

Insert the syringe properly into the TEOLINK

To firmly engage the syringe, pull the rubber stopper down until you hear a «click»


3 TEOLINK connection

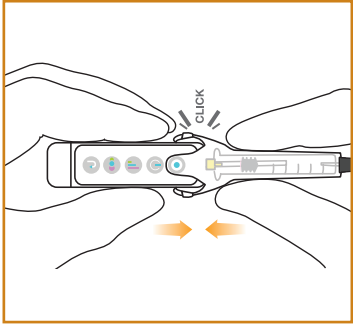


Engage the TEOLINK in the TEOSYAL®PEN

The «click» of the anchorage fins indicates that TEOSYAL®PEN is well engaged

4 Setting up the needle/cannula

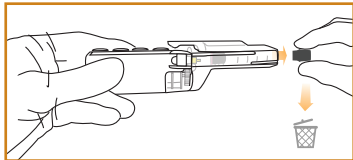


Remove and discard the cap from the syringe

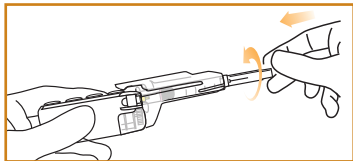


Place the needle or the cannula to the syringe

Remove the cap from the needle or the cannula and retain it for dismantling

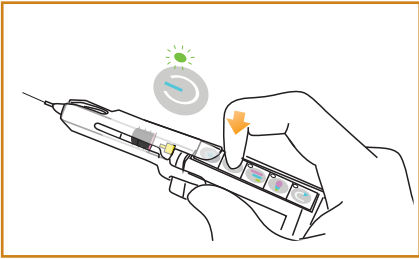


Turn on the device by pressing the ON/OFF button

The green light indicates that the device is switched on

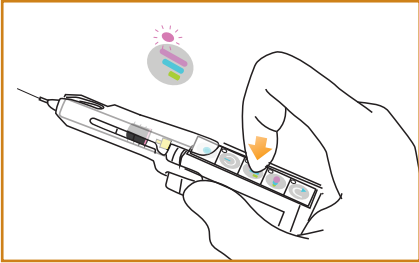

Select the «flow mode» and set-up the fast delivery speed

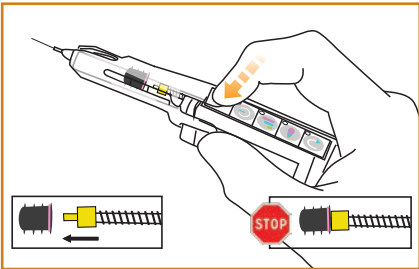

Move forward the piston by pressing and holding the activation button until reaching the rubber

The piston is correctly set when the head (in yellow) touches the rubber, then you have to stop pressing the activation button

Select the low speed «flow mode» to bring the gel until the end of the needle or cannula

According to preferences, two modes of injection exist:

Injection «flow mode»

The «flow mode» enables all injection techniques

Choose the speed according to the objectives

- Slow speed for superficial wrinkles and delicate zones
- Medium speed for superficial to moderately deep wrinkles
- Fast speed for deep injections

With this mode, it is also possible to realize multi-bolus injections by pressing alternatively the activation button

Injection «drop mode»

Select drop size according to the objectives

- Small drop
- Medium drop
- Large drop

At the start of slow injection speed, there is a half a second delay in the action

Use medium speed while becoming familiar with the pen before advancing to fast speed

During injection with «drop mode», press and hold the activation button, TEOSYAL®PEN automatically releases the desired drop size

Wait a few seconds after stopping TEOSYAL®PEN before removing the needle/the cannula from the skin

Do not forget that the flow also depends on the product fluidity, the selected injection speed and resistance from tissues

After the injection, press the rewind button to bring back the rod and turn off the device

TEOSYAL®PEN switches off automatically after approximately 4 minutes of inaction